


Datasheet

HBPR-450M-10K-IN-FST

High-Speed Balanced Photoreceiver

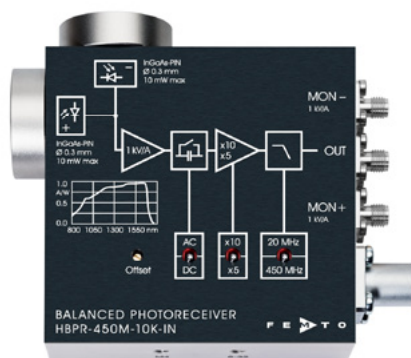


Features

- Bandwidth DC to 450 MHz
- Common-Mode Rejection Ratio (CMRR) 35 dB typ.
- InGaAs-PIN detectors, 0.3 mm active diameter
- Spectral range 800 – 1700 nm
- Very low NEP, down to 6.5 pW/√Hz
- Transimpedance gain switchable 5×10^3 V/A, 10×10^3 V/A
- High dynamic input range up to 2×10 mW balanced optical power
- Fast monitor outputs with 10 MHz bandwidth and 1×10^3 V/A gain
- Switchable low pass filter for minimizing wideband noise
- Free-space input 1.035"-40 threaded
- UNC 8-32 and M4 tapped holes for mounting on standard posts with metric and imperial thread

Applications

- Spectroscopy
- Heterodyne detection
- Optical coherence tomography (OCT)
- Optical delay measurement
- Differential optical front-end for oscilloscopes, spectrum analyzers, A/D converters and RF lock-in amplifiers

Block Diagram

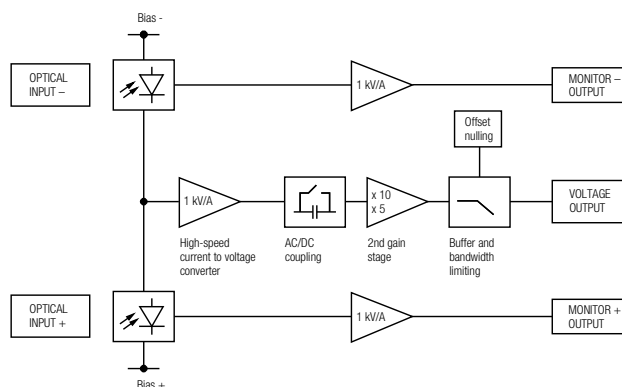

BS-HBPR_R2

SOPHISTICATED TOOLS FOR SIGNAL RECOVERY

FEMTO

HBPR-450M-10K-IN-FST_R3/TH_JMa/05MAR2024

Datasheet

HBPR-450M-10K-IN-FST

High-Speed Balanced Photoreceiver

Intended Use	The HBPR-450M-10K-IN-FST photoreceiver consists of a combination of two anti-parallel connected photodiodes with a subsequent low-noise transimpedance amplifier. It is designed for fast conversion of the tiny difference of two optical signals into an equivalent output voltage. Operation is mostly self-explanatory. If in doubt, consult this document or contact support@femto.de. For safe operation, please refer to the damage thresholds specified in the "Absolute Maximum Ratings", "Temperature Range" and "Power Supply" sections of this document. The operating environment must be free of smoke, dust, grease, oil, condensing moisture, and other contaminants that could affect the operation or performance.												
Application Notes	The damage threshold of 12 mW for each photodiode mentioned in the "Absolute Maximum Ratings" section applies to reasonably homogeneous illumination of the photodiodes. Extreme focusing of the light beam can lead to damage to the photodiodes, even at significantly lower light power. To achieve optimum performance, it is recommended that the CW light intensity at both inputs be well balanced. The monitor outputs can be used for continuous balance control. For setups with arbitrarily varying CW offset, the photoreceiver's AC mode can be helpful. Using AC mode increases the CW offset range to 1.3 mW (@ 1550 nm), regardless of the gain setting.												
Available Version	HBPR-450M-10K-IN-FST  1.035"-40 threaded flanges with internally threaded coupler rings mounted (outer dia. 30 mm), for free space applications, compatible with many optical standard accessories												
Related Models	Various free space or fiber coupled HBPR models, with bandwidth up to 500 MHz, in the spectral range from 320 nm to 1700 nm are available.												
Si Versions	Fiber-coupled with fix/permanent FC fiber connectors <table> <tr> <td>HBPR-100M-60K-SI-FC</td><td>Si-PIN Ø 0.8 mm, DC – 100 MHz, 320 – 1000 nm, CMRR 50 dB, gain $2.0 \times 10^4 / 6.0 \times 10^4$ V/A switchable</td></tr> <tr> <td>HBPR-200M-30K-SI-FC</td><td>Si-PIN Ø 0.8 mm, DC – 200 MHz, 320 – 1000 nm, CMRR 45 dB, gain $1.0 \times 10^4 / 3.0 \times 10^4$ V/A switchable</td></tr> <tr> <td>HBPR-500M-10K-SI-FC</td><td>Si-PIN Ø 0.4 mm, DC – 500 MHz, 320 – 1000 nm, CMRR 40 dB, gain $5.0 \times 10^3 / 10.0 \times 10^3$ V/A switchable</td></tr> </table> Free space versions with 1.035"-40 threaded flanges <table> <tr> <td>HBPR-100M-60K-SI-FST</td><td>Si-PIN Ø 0.8 mm, DC – 100 MHz, 320 – 1000 nm, CMRR 50 dB, gain $2.0 \times 10^4 / 6.0 \times 10^4$ V/A switchable</td></tr> <tr> <td>HBPR-200M-30K-SI-FST</td><td>Si-PIN Ø 0.8 mm, DC – 200 MHz, 320 – 1000 nm, CMRR 45 dB, gain $1.0 \times 10^4 / 3.0 \times 10^4$ V/A switchable</td></tr> <tr> <td>HBPR-500M-10K-SI-FST</td><td>Si-PIN Ø 0.4 mm, DC – 500 MHz, 320 – 1000 nm, CMRR 40 dB, gain $5.0 \times 10^3 / 10.0 \times 10^3$ V/A switchable</td></tr> </table>	HBPR-100M-60K-SI-FC	Si-PIN Ø 0.8 mm, DC – 100 MHz, 320 – 1000 nm, CMRR 50 dB, gain $2.0 \times 10^4 / 6.0 \times 10^4$ V/A switchable	HBPR-200M-30K-SI-FC	Si-PIN Ø 0.8 mm, DC – 200 MHz, 320 – 1000 nm, CMRR 45 dB, gain $1.0 \times 10^4 / 3.0 \times 10^4$ V/A switchable	HBPR-500M-10K-SI-FC	Si-PIN Ø 0.4 mm, DC – 500 MHz, 320 – 1000 nm, CMRR 40 dB, gain $5.0 \times 10^3 / 10.0 \times 10^3$ V/A switchable	HBPR-100M-60K-SI-FST	Si-PIN Ø 0.8 mm, DC – 100 MHz, 320 – 1000 nm, CMRR 50 dB, gain $2.0 \times 10^4 / 6.0 \times 10^4$ V/A switchable	HBPR-200M-30K-SI-FST	Si-PIN Ø 0.8 mm, DC – 200 MHz, 320 – 1000 nm, CMRR 45 dB, gain $1.0 \times 10^4 / 3.0 \times 10^4$ V/A switchable	HBPR-500M-10K-SI-FST	Si-PIN Ø 0.4 mm, DC – 500 MHz, 320 – 1000 nm, CMRR 40 dB, gain $5.0 \times 10^3 / 10.0 \times 10^3$ V/A switchable
HBPR-100M-60K-SI-FC	Si-PIN Ø 0.8 mm, DC – 100 MHz, 320 – 1000 nm, CMRR 50 dB, gain $2.0 \times 10^4 / 6.0 \times 10^4$ V/A switchable												
HBPR-200M-30K-SI-FC	Si-PIN Ø 0.8 mm, DC – 200 MHz, 320 – 1000 nm, CMRR 45 dB, gain $1.0 \times 10^4 / 3.0 \times 10^4$ V/A switchable												
HBPR-500M-10K-SI-FC	Si-PIN Ø 0.4 mm, DC – 500 MHz, 320 – 1000 nm, CMRR 40 dB, gain $5.0 \times 10^3 / 10.0 \times 10^3$ V/A switchable												
HBPR-100M-60K-SI-FST	Si-PIN Ø 0.8 mm, DC – 100 MHz, 320 – 1000 nm, CMRR 50 dB, gain $2.0 \times 10^4 / 6.0 \times 10^4$ V/A switchable												
HBPR-200M-30K-SI-FST	Si-PIN Ø 0.8 mm, DC – 200 MHz, 320 – 1000 nm, CMRR 45 dB, gain $1.0 \times 10^4 / 3.0 \times 10^4$ V/A switchable												
HBPR-500M-10K-SI-FST	Si-PIN Ø 0.4 mm, DC – 500 MHz, 320 – 1000 nm, CMRR 40 dB, gain $5.0 \times 10^3 / 10.0 \times 10^3$ V/A switchable												

SOPHISTICATED TOOLS FOR SIGNAL RECOVERY

F E M T O

HBPR-450M-10K-IN-FST_R3/TH_JMa/05MAR2024

Datasheet

HBPR-450M-10K-IN-FST

High-Speed Balanced Photoreceiver

Related Models (continued)

InGaAs Versions

Fiber-coupled with fix/permanent FC fiber connectors (ball lense coupled)

HBPR-100M-60K-IN-FC InGaAs-PIN $\varnothing$ 0.08 mm, DC – 100 MHz, 900 – 1700 nm, CMRR 55 dB, gain $2.0 \times 10^4 / 6.0 \times 10^4$ V/A switchable

HBPR-200M-30K-IN-FC InGaAs-PIN $\varnothing$ 0.08 mm, DC – 200 MHz, 900 – 1700 nm, CMRR 50 dB, gain $1.0 \times 10^4 / 3.0 \times 10^4$ V/A switchable

HBPR-500M-10K-IN-FC InGaAs-PIN $\varnothing$ 0.08 mm, DC – 500 MHz, 900 – 1700 nm, CMRR 45 dB, gain $5.0 \times 10^3 / 10.0 \times 10^3$ V/A switchable

Free space versions with 1.035"-40 threaded flanges

HBPR-100M-60K-IN-FST InGaAs-PIN $\varnothing$ 0.3 mm, DC – 100 MHz, 800 – 1700 nm, CMRR 50 dB, gain $2.0 \times 10^4 / 6.0 \times 10^4$ V/A switchable

HBPR-200M-30K-IN-FST InGaAs-PIN $\varnothing$ 0.3 mm, DC – 200 MHz, 800 – 1700 nm, CMRR 45 dB, gain $1.0 \times 10^4 / 3.0 \times 10^4$ V/A switchable

Available Accessory

PS-15-25-L



Power Supply
Input: 100 – 240 VAC
Output: ± 15 VDC

Specifications

Test conditions

$V_s = \pm 15$ V, $T_A = 25$ °C, output load impedance 50 Ω , warm-up 20 minutes (min. 10 minutes recommended), monitor outputs terminated with 1 M Ω

Gain

Transimpedance gain

5×10^3 V/A (@ 2nd gain $\times 5$, 50 Ω load)
 10×10^3 V/A (@ 2nd gain $\times 10$, 50 Ω load)

Gain accuracy

± 1 % electrical

Conversion gain

4.75×10^3 V/W typ. (@ 2nd gain $\times 5$, 1550 nm, 50 Ω load)
 9.50×10^3 V/W typ. (@ 2nd gain $\times 10$, 1550 nm, 50 Ω load)

Common mode rejection ratio (CMRR)

50 dB typ. ($f \leq 100$ MHz)
35 dB typ. ($f \leq 450$ MHz)

Frequency Response

Lower cut-off frequency

DC / 10 Hz, switchable

Upper cut-off frequency (–3 dB)

450 MHz / 20 MHz, switchable

Time Response

Rise/fall time (10 % – 90 %)

0.95 ns (@ 2nd gain $\times 5$); 1.02 ns (@ 2nd gain $\times 10$)
17.5 ns (@ bandwidth set to 20 MHz)

Input

Noise equivalent power (NEP)

minimum 6.5 pW/√Hz (@ 1550 nm)
6.9 pW/√Hz (@ 1550 nm, 20 MHz)
20.0 pW/√Hz (@ 1550 nm, 200 MHz)
38.0 pW/√Hz (@ 1550 nm, 450 MHz)

Maximum differential CW power (for linear amplification)

210 μ W (@ 2nd gain $\times 5$, DC-coupled, 1550 nm)
105 μ W (@ 2nd gain $\times 10$, DC-coupled, 1550 nm)

Max. optical CW balanced power (common mode power)

1.3 mW (@ AC-coupled, 1550 nm)
10 mW (on each photodiode, @ 1550 nm)

Monitor optical saturation power (limited by linear amplification)

10.5 mW (@ 1550 nm)

Detector

Detector type

InGaAs-PIN photodiode

Active area

$\varnothing$ 300 μ m

Spectral range

800 – 1700 nm

Sensitivity

0.95 A/W typ. (@ 1550 nm)

SOPHISTICATED TOOLS FOR SIGNAL RECOVERY

F E M T O


HBPR-450M-10K-IN-FST_R3/TH_JMa/05MAR2024

Datasheet
HBPR-450M-10K-IN-FST
High-Speed Balanced Photoreceiver
Specifications (continued)

Output	Output voltage range Max. output voltage Offset voltage compensation Output impedance Slew rate Max. output current Output reflection S22 Output noise (typ.)	± 1.0 V (@ 50 Ω load) for linear operation and low harmonic distortion ± 2.0 V (@ 50 Ω load) ± 100 mV typ., adjustable by offset potentiometer 50 Ω (terminate with 50 Ω load) 2800 V/ μ s 70 mA -30 dB @ < 100 MHz -20 dB @ < 800 MHz 2.4 mV RMS (16 mV peak-peak) (@ 2 nd gain $\times 5$) 4.6 mV RMS (31 mV peak-peak) (@ 2 nd gain $\times 10$) 0.25 mV RMS (1.7 mV peak-peak) (@ 2 nd gain $\times 5$, BW 20 MHz) 0.4 mV RMS (2.7 mV peak-peak) (@ 2 nd gain $\times 10$, BW 20 MHz) (@ 50 Ω load, no signal on detectors, measurement bandwidth 2 GHz)
Monitor Outputs	Gain Voltage range Output impedance Max. output current Bandwidth Output noise	1×10^3 V/A (@ ≥ 100 k Ω load) 0 ... +10 V (@ ≥ 100 k Ω load) 50 Ω (terminate with ≥ 100 k Ω load) 30 mA typ. DC – 10 MHz 0.6 mV RMS (4 mV peak-peak) (@ 100 k Ω load, no signal on detectors, measurement bandwidth 200 MHz)
Power Supply	Supply voltage Supply current	± 15 V (± 14.5 V ... ± 16.5 V) -90 / $+120$ mA typ. (depends on operating conditions, recommended power supply capability min. ± 200 mA)
Optical Input Connector	Material FST flange Material FST coupler ring	1.4305 stainless steel, nickel-plated 1.4305 stainless steel, glass bead blasted
Case	Weight Material	410 g (0.9 lbs) including coupler rings AlMg3Mn, nickel-plated
Temperature Range	Storage temperature Operating temperature	-40 °C ... $+85$ °C 0 °C ... $+60$ °C

Absolute Maximum Ratings	Optical input power (CW) Power supply voltage	12 mW (on each photodiode) ± 20 V
--------------------------	--	--

Connectors	Inputs Outputs Power supply	1.035"-40 threaded flanges for free space applications and for use with various types of optical standard accessories SMA jacks (female) LEMO® series 1S, 3-pin fixed socket (mating plug type: FFA.1S.303.CLAC52)
------------	-----------------------------------	--

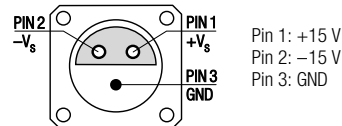

Scope of Delivery	HBPR-450M-10K-IN-FST, 2 $\times$ threaded coupler ring, Lemo® 3-pin connector, 3 $\times$ adapter SMA (male) to BNC (female), datasheet	
-------------------	--	--

SOPHISTICATED TOOLS FOR SIGNAL RECOVERY
F E M T O

HBPR-450M-10K-IN-FST_R3/TH_JMa/05MAR2024


Datasheet
HBPR-450M-10K-IN-FST
High-Speed Balanced Photoreceiver
Ordering Information

HBPR-450M-10K-IN-FST

1.035"-40 threaded flanges for free space applications and for use with various types of optical standard accessories

Spectral Response

HBPR-450M-10K-IN-FST

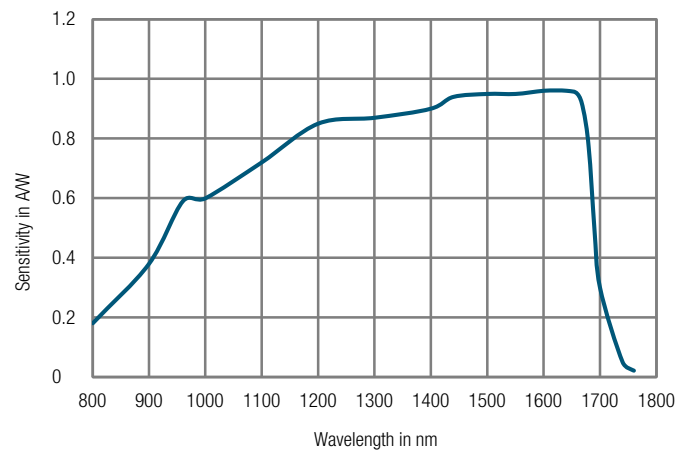

DB-Sens-HBPR-IN-FST_R2

SOPHISTICATED TOOLS FOR SIGNAL RECOVERY
F E M T O

HBPR-450M-10K-IN-FST_R3/TH_JMa/05MAR2024


Datasheet

HBPR-450M-10K-IN-FST

High-Speed Balanced Photoreceiver

Dimensions

HBPR-450M-10K-IN-FST

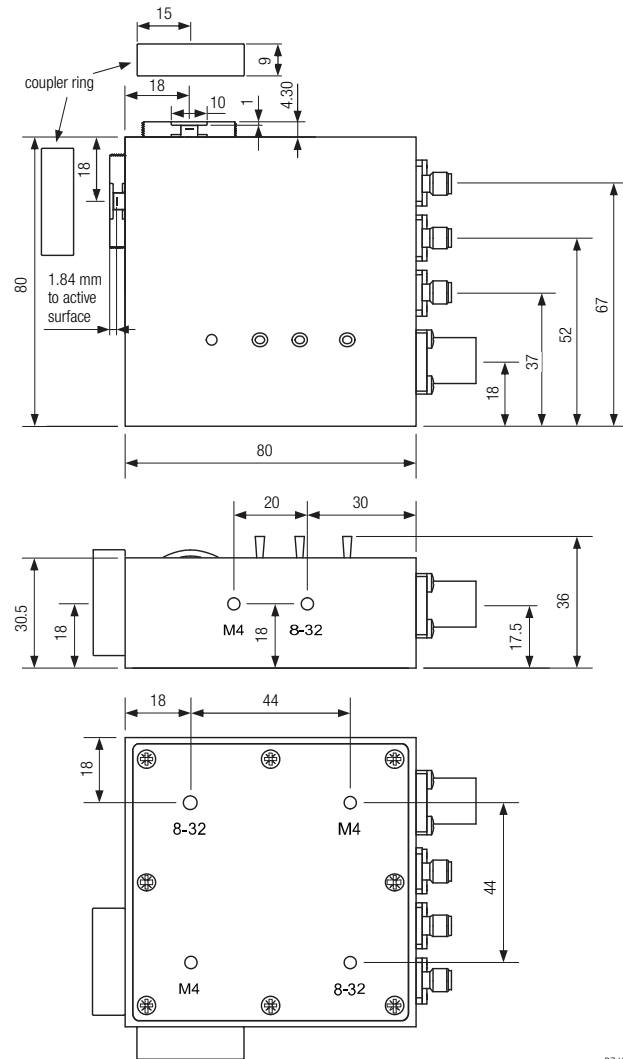

all dimensions in mm unless otherwise noted

DZ-HBPR_FST_R2

The bottom plate may be rotated to match the appropriate mounting thread to the optical axis by unscrewing the 8 screws.

Specifications are subject to change without notice. Information provided herein is believed to be accurate and reliable. However, no responsibility is assumed by FEMTO Messtechnik GmbH for its use, nor for any infringement of patents or other rights of third parties which may result from its use. No license is granted by implication or otherwise under any patent or patent rights of FEMTO Messtechnik GmbH. Product names mentioned may also be trademarks used here for identification purposes only.

© by FEMTO Messtechnik GmbH · Printed in Germany

SOPHISTICATED TOOLS FOR SIGNAL RECOVERY

F E M T O

HBPR-450M-10K-IN-FST_R3/TH_JMa/05MAR2024